

**PCT**WORLD INTELLECTUAL PROPERTY ORGANIZATION  
International Bureau

## INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

<b>(51) International Patent Classification <sup>6</sup> :</b> <b>G01N 33/80, 33/543</b>	<b>A1</b>	<b>(11) International Publication Number:</b> <b>WO 95/31731</b> <b>(43) International Publication Date:</b> 23 November 1995 (23.11.95)
<b>(21) International Application Number:</b> PCT/US95/05982 <b>(22) International Filing Date:</b> 15 May 1995 (15.05.95)  <b>(30) Priority Data:</b> 08/243,296 17 May 1994 (17.05.94) US  <b>(71) Applicant:</b> GAMMA BIOLOGICALS, INC. [US/US]; 3700 Mangum Road, Houston, TX 77092-5497 (US). <b>(72) Inventors:</b> FRAMB, Thomas, H.; 8915 Herts Road, Spring, TX 77379 (US). HATCHER, David, E.; 15902 Champion Drive, Houston, TX 77079 (US). MOULDS, John, J.; 10235 Oak Point Drive, Houston, TX 77043 (US). <b>(74) Agent:</b> EDWARDS, Stephen, E.; Arnold, White & Durkee, P.O. Box 4433, Houston, TX 77210 (US).		<b>(81) Designated States:</b> AM, AT, AU, BB, BG, BR, BY, CA, CH, CN, CZ, DE, DK, EE, ES, FI, GB, GE, HU, IS, JP, KE, KG, KP, KR, KZ, LK, LR, LT, LU, LV, MD, MG, MN, MW, MX, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, TJ, TT, UA, UG, UZ, VN, European patent (AT, BE, CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, ML, MR, NE, SN, TD, TG), ARIPO patent (KE, MW, SD, SZ, UG).  <b>Published</b> <i>With international search report.</i> <i>Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.</i>
<b>(54) Title:</b> METHOD AND APPARATUS USEFUL FOR DETECTING BLOODGROUP ANTIGENS AND ANTIBODIES  <b>(57) Abstract</b>  The direct assay comprises adding a sample of erythrocytes to a reaction tube charged with a column of immunoreactive particles having an immunoglobulin binding ligand selected from the group consisting of Protein A, Protein G, Protein A/G or a universal kappa light chain binding protein coupled to the surface of the particles. Antibodies specific for bloodgroup antigens tested for are coupled to the ligand on the particles. The indirect assay comprises obtaining either a sample of erythrocytes or a sample of blood serum to be tested and mixing the erythrocytes or serum with a known antibody or antigen reagent, depending on whether antigens or antibodies are being tested for. The mixture is incubated in a reaction tube above a column of immunoreactive particles having immunoglobulin binding ligands selected from the group consisting of Protein A, Protein G, Protein A/G or a universal kappa light chain binding protein coupled to the surface of the particles. A positive response is indicated by binding of the erythrocytes to the immunoreactive particles in the upper one-half of the reaction tube. A negative response is indicated by the erythrocytes being forced to the bottom of the reaction tube in a pellet.		

BEST AVAILABLE COPY

**FOR THE PURPOSES OF INFORMATION ONLY**

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	GB	United Kingdom	MR	Mauritania
AU	Australia	GE	Georgia	MW	Malawi
BB	Barbados	GN	Guinea	NE	Niger
BE	Belgium	GR	Greece	NL	Netherlands
BF	Burkina Faso	HU	Hungary	NO	Norway
BG	Bulgaria	IE	Ireland	NZ	New Zealand
BJ	Benin	IT	Italy	PL	Poland
BR	Brazil	JP	Japan	PT	Portugal
BY	Belarus	KE	Kenya	RO	Romania
CA	Canada	KG	Kyrgyzstan	RU	Russian Federation
CF	Central African Republic	KP	Democratic People's Republic of Korea	SD	Sudan
CG	Congo	KR	Republic of Korea	SE	Sweden
CH	Switzerland	KZ	Kazakhstan	SI	Slovenia
CI	Côte d'Ivoire	LI	Liechtenstein	SK	Slovakia
CM	Cameroon	LK	Sri Lanka	SN	Senegal
CN	China	LU	Luxembourg	TD	Chad
CS	Czechoslovakia	LV	Latvia	TG	Togo
CZ	Czech Republic	MC	Monaco	TJ	Tajikistan
DE	Germany	MD	Republic of Moldova	TT	Trinidad and Tobago
DK	Denmark	MG	Madagascar	UA	Ukraine
ES	Spain	ML	Mali	US	United States of America
FI	Finland	MN	Mongolia	UZ	Uzbekistan
FR	France			VN	Viet Nam
GA	Gabon				

5

10                   **METHOD AND APPARATUS USEFUL FOR DETECTING**  
                      **BLOODGROUP ANTIGENS AND ANTIBODIES**

                  The present invention relates generally to the field of immunoassay testing  
and more particularly to the detection of bloodgroup antigens and antibodies.

15

                  The testing of blood for various antigens and antibodies has become  
commonplace. For example, before a patient receives a blood transfusion, a  
number of laboratory tests must be performed on a sample of both the patient's  
and the donor's blood. Specifically, the blood must be typed, to determine its  
20    ABO and Rh D bloodgroup. Before a transfusion can proceed, the patient's serum  
must be tested for compatibility, or cross-matched, to determine whether the  
serum contains antibodies to antigens present on the erythrocytes (red blood cells)  
of a specific donor.

25

                  Additionally, there are an increasing number of antibodies that have been  
determined to be present in blood serum. Determining the presence or absence of  
certain antibodies can have great importance, both in the diagnosis and the  
treatment of certain disorders.

30

                  Blood typing, compatibility testing, and antibody screening are usually done  
through some type of agglutination immunoassay. Agglutination assays, due to  
their simplicity and wide range of applications, have become one of the most  
common methods used in blood testing. Agglutination assays can be performed  
relatively simply and can be detected or read visually without resort to expensive

detection equipment. In a classical simple direct agglutination assay, e.g., agglutination of group A erythrocytes by anti-A sera, cells bearing a group A antigen are agglutinated directly by antibody.

5           One serious drawback to this type of testing is the tendency for erythrocytes possessing weak expression of a particular bloodgroup antigen being tested for to not agglutinate or clump to an extent so that proper visual determination of the positive response is possible. Moreover, in some cases, IgG class antibodies bind to the erythrocytes with high affinity but fail to induce  
10           agglutination. Addition of anti-IgG antibodies cross links the antibodies located on the erythrocytes and causes agglutination. An example of this type of reaction is the classical indirect Coombs test in which a test serum is first incubated with test erythrocytes, the cells are washed to remove excess unbound antibodies, and then mixed with an antiglobulin serum (Coombs reagent).

15           Lapierre et al. reported a variation of this agglutination assay that uses inert particles in the presence of a reagent that can be either an antibody or an antigen, for typing blood components that can be either erythrocytes or serum. *Transfusion* 30:109-113 (1990). For an indirect test, the reagent (e.g., serum) and blood  
20           components (e.g., erythrocytes) are mixed and incubated for a specified period of time. The mixture is centrifuged for about 10 minutes at 70 x g into a bed of inert particles (i.e., a dextran gel) in a medium containing an agglutinating reagent in solution, generally an anti-human IgG (Coombs reagent). The erythrocytes that agglutinate are trapped within the matrix of inert particles, indicating a positive  
25           reaction. In a negative reaction, the erythrocytes circumvent the particles and pellet at the bottom of the reaction tube.

          Because the Lapierre method indicates a positive response by trapping the agglutinated erythrocytes on the matrix of inert particles, weak positive reactions  
30           vary significantly in appearance. Generally, in a weak positive reaction, some of the erythrocytes will pellet in the bottom of the reaction tube and others will be

trapped in multiple clumps. The size of the clumps of agglutinated erythrocytes will be dependent upon the strength of the antibody/antigen reaction. Thus, where the reaction is weak, the clumps will be smaller and more likely to pass through the particle matrix. Although the Lapierre method offers an easy-to-perform test,  
5 weak positive reactions may still be difficult to read.

Accordingly, there remains room in this field for improved techniques in blood typing assays, in making an assay that is easier, more flexible and more predictable to read where there is a weak positive response. The present invention  
10 provides an improved, flexible, rapid, and accurate bloodgrouping system. The result of such a test is fewer false negative responses.

The method of the present invention is directed toward the detection of bloodgroup antigens and antibodies. It employs immunoreactive affinity  
15 chromatography techniques to detect bloodgroup antigens and antibodies.

This method uses a particle to which an antibody binding ligand is covalently bound. Thus, the particle is immunologically reactive ("immunoreactive"). The immunoreactive particle may be of various materials,  
20 although agarose is preferred. The ligand may be one of several immunoglobulin binding proteins, for example, Protein A, Protein G, Protein A/G or KappaLock™, although Protein G is preferred for applications in which antibody of the IgG class is used (e.g., in a direct test) or tested for (e.g., in an indirect test).

25 In direct assays, erythrocytes are layered over a bed of particles in a reaction tube. The immunoreactive particles have antibody specific for the antigen of interest bound to their surface through the ligand. The reaction tube is centrifuged and non-binding erythrocytes will pellet at the bottom of the tube. A positive reaction will cause the erythrocytes to bind to the particles in the upper  
30 portion of the reaction tube and will leave a distinct line of bound red cells in the upper one-half of the reaction tube. A negative response leaves substantially all of

the erythrocytes pelleted at the bottom of the reaction tube. In the case of a weak positive reaction, the erythrocytes with adequate quantity or quality of the tested for antigen will bind to the particles in the upper one-half of the reaction tube leaving a distinct zone. Most of the remainder of the erythrocytes will pellet at  
5 the bottom of the reaction tube. The binding of the erythrocytes to the particles is due to the affinity of the antigens on the red cells for the antibodies bound to the ligand on the particles.

In indirect assays, erythrocytes and serum are incubated together for a  
10 sufficient time for the antigens on the surface of the erythrocytes to react with the antibodies present in the serum. The erythrocyte/serum mixture is centrifuged, exposing the erythrocyte-bound antibodies to the immunoreactive particles. The centrifuging step is for a time sufficient to force any non-binding erythrocytes through the matrix of particles. A positive reaction leaves a distinct line of  
15 erythrocytes bound to the particles due to the affinity of the ligand for the antibodies that have become attached to the antigens on the erythrocytes. A negative response leaves the erythrocytes pelleted at the bottom of the reaction tube. In the case of a weak positive response, some of the erythrocytes will form a line in the upper one-half of the reaction tube and the remainder will pellet at the  
20 bottom of the reaction tube.

The reaction that binds the erythrocytes to the particles in the direct assay, i.e., an antigen/antibody reaction, is the same reaction that causes erythrocytes to agglutinate. Here, however, instead of merely forming clumps too large to pass  
25 through the spaces between the particles in the reaction tube, the cells bind directly to the particle and not to one another. Hence, the reaction of the present invention is an affinity reaction and not an agglutination reaction.

Similarly, in indirect assays, the erythrocytes become bound to the particles  
30 through the affinity that the ligand has for the antibody which has become affixed

to the erythrocyte antigen during the incubation period. There is no agglutination of the erythrocytes but rather an adherence of the erythrocytes to the particles.

Thus, the method of the present invention relies upon the affinity of the  
5 erythrocytes for the particle to indicate a positive response. Consequently, positive reactions are more pronounced and definite in the method of the present invention than in those of the prior art.

FIG. 1A is a side view of a strong positive reaction wherein the binding  
10 erythrocytes form a very distinct and broad line in the upper one-half of the reaction tube.

FIG. 1B is a side view of a weak positive reaction wherein an adequate  
quantity or quality of antigen/antibody reaction causes the erythrocytes to bind to  
15 the particles forming a distinct line in the upper one-half of the reaction tube. The erythrocytes not creating the quantity or quality of antigen/antibody reaction to cause the same binding will be forced to pellet in the bottom of the reaction tube.

FIG. 1C is a side view of a negative reaction wherein none of the  
20 erythrocytes create the antigen/antibody reaction necessary to bind the erythrocytes to the particles and as they did not bind to the particles are forced to pellet in the bottom of the reaction tube.

FIG. 2A is a perspective view of a preferred single unit 6-tube array.  
25

FIG. 2B is a top view of the same single unit 6-tube array.

FIG. 2C is a side view of the same single unit 6-tube array.

30 FIG. 2D is an end view of the same single unit 6-tube array showing a squared-off base.

FIG. 2E is an end view of the end opposite that shown in FIG. 2D and showing a slightly narrower base than seen in FIG. 2D.

The method of the present invention has two preferred embodiments. The  
5 first is a direct assay, in which an antibody against an erythrocyte antigen is coupled to a ligand that is covalently coupled to a particle within the reaction tube. The erythrocyte is then bound to the particles through the ligand-bound antibody. The second embodiment is an indirect assay, where erythrocytes are exposed to serum antibodies and the mixture is incubated to bind the antibodies to the  
10 erythrocytes, and the erythrocyte-antibody complexes are then bound to the ligand on the particles.

The particles used in the method of the present invention may be any particle used in agglutination assays to which a ligand may be coupled.  
15 Preferably, however, the particles are of agarose, Sepharose or latex, although other types of particles to which ligand may be coupled are also included within the scope of the invention. These particles are generally in the form of beads. Most preferred are agarose and Sepharose. These particles are in the form of beaded gels.

20 A number of ligands are known that bind immunoglobulin molecules and may be covalently coupled to the particles, for example Protein A, Protein G, Protein A/G and KappaLock™. Protein G is a particularly preferred ligand for use in assays where IgG immunoglobulins are used or tested for. One reason that  
25 Protein G is preferred is that it has a greater affinity than Protein A for most IgG immunoglobulins. Protein G also binds with a significantly greater affinity than Protein A to certain subclasses of IgG, e.g., human IgG3, mouse IgG1 and rat IgG2a. Protein G does not bind to human IgM, IgA and IgD.

30 Protein G is a bacterial cell wall protein isolated and purified from group G streptococci. Protein G binds to mammalian IgG immunoglobulins through their



Fc portion. Because Protein G only binds the Fc portion of IgG immunoglobulins, the antibody portion of the immunoglobulin remains available for reaction with its corresponding antigen, yet the immunoglobulin remains bound to the particle.

Native Protein G has been sequenced by DNA analysis. From the DNA analysis,  
5 two IgG binding domains and sites for albumin and cell surface binding have been identified.

ImmunoPure® Immobilized Protein G is a commercially available particle product having Protein G immobilized on the surface of agarose gel bead particles.  
10 This product is available from Pierce of Rockford, Illinois. The immobilized Protein G has been genetically engineered to remove the albumin and cell surface binding regions and thereby minimize the binding of anything other than immunoglobulins.

15 ImmunoPure® Immobilized Protein G consists of a recombinant Protein G covalently linked (glutaraldehyde activation of the beads) to cross-linked 6% beaded agarose. The material is supplied in a 50% slurry. The material can bind 11 mg of human IgG per ml of gel.

20 Protein A is a cell wall component produced by several strains of *Staphylococcus aureus*. Protein A is capable of specifically binding the Fc region of immunoglobulin molecules, especially IgG. The Protein A molecule has four high affinity binding sites that are capable of interacting with the Fc region from IgG of several species. Protein A interacts with some IgG subgroups and not with  
25 others. For example, human IgG1, IgG2 and IgG4 bind strongly while IgG3 does not bind. And, there are also some instances in which monoclonal antibodies do not bind to Protein A.

Immobilized Protein A is also available commercially from Pierce. This  
30 immobilized Protein A is a highly purified Protein A, covalently coupled to

cross-linked beaded agarose. The typical binding capacity for this immobilized Protein A is 12-15 mg of human IgG per milliliter of gel.

Protein A/G is a genetically engineered protein that combines the IgG  
5 binding profiles of both Protein A and Protein G. Protein A/G is a gene fusion product secreted from a non-pathogenic form of *Bacillus*. This genetically engineered Protein A/G is designed to contain four Fc binding domains from Protein A and two from Protein G.

10 Protein A/G binds to all human IgG subclasses. In addition, it binds to IgA, IgE, IgM and to IgD but to a lesser extent to IgD. Thus, Protein A/G may be a preferred ligand in tests for or tests using non-IgG class immunoglobulins.

Pierce also offers an immobilized Protein A/G covalently coupled to beaded  
15 agarose under the trade name ImmunoPure® Immobilized Protein A/G.

KappaLock™ is a universal kappa light chain binding protein available from Aaston, Inc., 12 Falmouth Road, Wellesley, Massachusetts. It has been  
20 genetically engineered from the DNA of a strain of *Peptostreptococcus*. This protein binds to the kappa region of the light chain of all antibody types. KappaLock™ has been genetically engineered to delete the albumin and cell wall binding regions of the naturally occurring bacterial protein. The resultant engineered protein has four antibody binding domains and specifically does not  
25 bind to heavy chains or to the Fc region of immunoglobulins. Because kappa light chains are shared among different classes of antibodies, KappaLock™ will bind to antibodies having a kappa light chains regardless of heavy chain class.

KappaLock™ may be immobilized on various supports, particularly on agarose beads. Immobilized KappaLock™ will capture mouse IgG, rabbit IgG,  
30 human IgG, human IgA and human IgM.

All of the preferred ligands may be covalently bound to a solid-phase matrix such as agarose beads (e.g., Sepharose Pharmacia) using known techniques, for example as described by Hearn et al., *Methods in Enzymology* Vol. 35:102-117 (1987). Generally, the beads are first activated by a chemical agent, such as  
5 glutaraldehyde, carbonyldiimidazole, cyanogen bromide hydroxy-succinimide, tosyl chloride or the like. The chosen ligand is then covalently attached to the beads, resulting in an extremely stable linkage of the ligand to the support.

The antibody or antigen used as a reagent is dependent upon the antibody  
10 or antigen that is being tested for. The number of blood antigens and thus, antibodies that have been identified is very large, with more antigens and antibodies continually being determined. The International Society of Blood Transfusion has published an article titled *Bloodgroup Terminology 1990*, Vox. Sang. 58:152-169 (1990), which is expressly incorporated herein by reference. A  
15 non-exclusive list of red cell antigens for which the present method would be suitable is found on page 153. However, the following antibodies and antigens are particularly preferred: A, B, D, C, c, C<sup>w</sup>, E, e, K, Fy<sup>a</sup>, Fy<sup>b</sup>, Jk<sup>a</sup>, Jk<sup>b</sup>, S and s.

In adapting the method for use in testing for the various antibodies and  
20 antigens of interest, one should select a ligand capable of binding to the isotype of the antibody that is used or tested for or, alternatively, one may use a bridging antibody, e.g., an IgG anti-IgM, for an IgM antibody. Thus, an IgG anti-IgM antibody would be coupled to the ligand as a "bridge" and an IgM antibody would bind to the IgG anti-IgM antibody.

25 For direct assay, the immunoreactive particles may be prepared as follows: the particles with covalently coupled ligand are first centrifuged and the supernate is removed and discarded. An appropriate amount of properly diluted antibody reagent is gently mixed with the particles, and the mixture is allowed to stand for  
30 about 15 minutes so as to allow the antibody to bind to the ligand on the particles. The amount of diluted antibody reagent added to the particles will depend upon the

source and concentration of the antibody to be used. The mixture is centrifuged, and the supernate is removed and replaced by a non-antibody containing solution, for example a buffered solution containing dextran. Thus, the immunoreactive particle suspension will generally contain little or no free antibody in solution.

5

An appropriate amount of the immunoreactive particle suspension, e.g., 20-150 $\mu$ l, is pipetted into the reaction tube, and the particles are allowed to settle to the bottom of the tube, resulting in a substantially particle free zone of dextran-buffer covering the particles. In some embodiments, a number of neutral particles, i.e., those not having ligand bound to the surface of the particles, may be mixed with the immunoreactive particles in order to make the assay less expensive to perform.

Generally, the erythrocytes to be tested will have been drawn from a patient and submitted for determination of the presence of certain bloodgroup antigens. The erythrocytes are separated from the whole blood by traditional techniques, and a suspension of erythrocytes in the range of about 0.5-1.0% (w/v) is prepared in a low ionic strength solution. Gamma N-HANCE<sup>®</sup> by Gamma Biologicals, Inc., is preferred for this solution.

20

An appropriate amount of a suspension of erythrocytes, e.g., 15-75 $\mu$ l of a 0.5-1% suspension, is added to the reaction tube and the reaction tube is centrifuged for a time sufficient to force the non-binding erythrocytes to the bottom. Preferably, the centrifuge has a rotor adapted so that the centrifugal force generated by the centrifuge acts along the long axis of the reaction tube. Consequently, the centrifugal force generated during the centrifuging step will operate to more forcefully urge the non-binding erythrocytes to pellet in the bottom of the reaction tube.

30

Those of ordinary skill in the art will recognize that the centrifugation conditions used will depend on a number of factors, e.g., the design and type of

-11-

the reaction vessel used and the centrifuge used. A two- or three-step centrifugation procedure is preferred. Generally, the initial centrifugation step is for a short period of time ranging from about 10-30 seconds at a high centrifugal force ranging from about 800-1200g in order to force the erythrocytes through the dextran solution. A second centrifugation step is used to bring the erythrocytes into contact with the immunoreactive particles. This step can be performed at a lower centrifugal force, ranging from about 300-600g for about 15-45 seconds. A third centrifugation step is used to force the unbound erythrocytes through the bead matrix. The third centrifugation may be performed at a high centrifugal force ranging from about 800-1200g for a period of about 45-90 seconds. The third centrifugation step considerably reduces the time needed to perform the separation of binding and non-binding erythrocytes. Separation is also possible by extending the duration of the second centrifugation step for about 5-10 minutes.

In the indirect assay, an appropriate amount of immunoreactive particles, e.g., 15-150 $\mu$ l, having the ligand bound to the surface of the particles is obtained and an appropriate amount of a diluent, such as buffered saline containing dextran, e.g., 30-300 $\mu$ l, is mixed with the particles. When performing an indirect assay for an antibody screening, it may be advantageous to use a dextran/saline buffer having a pH of about 7. An appropriate amount of the mixture, e.g., 20-150 $\mu$ l, is pipetted into a reaction tube and the particles are allowed to settle for approximately 15 minutes, resulting in a substantially particle free zone of dextran-buffer covering the particles. In some embodiments, a number of neutral particles, i.e., those not having ligand bound to the surface of the particles, may be mixed with the immunoreactive particles in order to make the assay less expensive to perform.

An appropriate amount of potentiator, e.g., 15-100 $\mu$ l, is added to the reaction tube. An appropriate amount of a suspension of erythrocytes, e.g., 15-75 $\mu$ l, is suspended in a suitable medium, such as a low ionic strength potentiator solution. The potentiator acts to enhance and increase the speed at

-12-

which the antigen or antibody binding occurs and creates a more profound positive reaction if one is to occur. A particularly preferred potentiator is a low ionic strength enhancement solution sold under the trade name Gamma N-HANCE® marketed by Gamma Biologicals, Inc., Houston, Texas. An appropriate amount  
5 of a suspension of erythrocytes to be tested, e.g., 15-75 $\mu$ l of a 0.5-1% suspension, is suspended in the potentiator solution, followed by the addition of an appropriate amount of test serum, e.g., 15-75 $\mu$ l. The reaction tube is then incubated under conditions sufficient for the antibody/antigen reaction to occur. Preferred conditions are incubation for approximately 5-30 minutes at about 37°C,  
10 and more preferred is an incubation time of about 10 minutes.

The reaction tube is then centrifuged for a time sufficient to force non-binding erythrocytes through the spaces between the particles, forcing them to pellet at the bottom of the reaction tube.

15 In comparing the indirect assay with the direct assay, it will be noted that in the indirect assay, the antibodies are added to the erythrocytes and that mixture incubated prior to the centrifuging step. Thus, the antibodies and antigens have the opportunity to react before the erythrocytes are exposed to the particles and the  
20 attached ligand. In the direct assay, the antibodies are coupled to the ligand on the particles and the erythrocytes are not exposed to the antibodies until the erythrocytes are forced through the spaces between particles by the centrifuge.

In the indirect assay, the erythrocytes that react with antibodies in the test  
25 serum bind to the particles through the interaction between the antibodies that attach to the antigens on the erythrocytes and the ligand on the particle. The detection of a positive reaction is the same as in the direct assay, in that a positive response results in a strong line of bound erythrocytes in the upper one-half of the reaction tube, a weak positive response results in a line of bound erythrocytes in  
30 the upper one-half of the reaction tube and generally some non-binding

erythrocytes pelleted at the bottom of the tube, and a negative response results in all the non-binding erythrocytes pelleted at the bottom of the reaction tube.

5 The reaction tube used in the method of the present invention is not critical. The most important factor to consider is that the tube provide an adequate column height of particles so that the affinity reaction can occur during centrifugation.

10 Microwells in microliter plates have been used successfully wherein a 3 to 4mm column of particles is present. Also useful are microfuge tubes and the like.

15 In a preferred embodiment, a plurality of reaction tubes can be joined to produce a single unit tube array, for example, as shown in FIG. 2A. FIG. 2A shows a single unit tube array comprising six reaction tubes 10. Each reaction tube is sealable at its top 12 so that the array may be sold prefilled with activated particles suitable for use in direct assays, indirect assays or a mixture thereof.

20 As seen in FIG. 2C, each reaction tube is substantially wider at the top, forming a reaction well 14 above the tube portion 16 of the reaction tube. The reaction well narrows in a tapering zone 18 so that during the centrifugation step, the erythrocytes present in the reaction well are directed into the tube portion of the reaction tube. When prefilled, the reaction tubes will be filled with immunoreactive particles and dextran solution into the tapering zone.

25 FIG. 2B shows that the tapering zone possesses four planar surfaces 20 with valleys 22 at the intersection of each of the planar surfaces. Each planar surface is broadest at the top of the zone where it meets the reaction well and then each surface gradually narrows through the zone until it reaches the reaction tube at the termination of the tapering zone.

30

Each reaction tube has a pointed terminus 24. This configuration causes non-binding erythrocytes to collect in the pointed terminus of the reaction tube in a pronounced pellet. The tube portion 16 of the reaction tube is to be dimensioned to provide a column of immunoreactive particles about 3-4mm in height. The array is designed to have a labelling surface 26 for affixing identifying information.

As seen in FIGS. 2D and 2E, each end of the single unit tube array is different, creating an index for orientation of these arrays during mass filling of the reaction tubes. Thus, one end of the array may be distinguished from the other.

The unit tube array may be manufactured using any suitable material, such as plastics and the like, but preferably is made of TPX®, a biologically and chemically inert plastic generically known as polymethylpentene. TPX® is supplied by Mitsui Petrochemicals (America) Ltd.

In one embodiment, e.g., for a direct test, the tubes in the array may be custom filled with a preselected variety of antibodies bound to the ligand. For example, a preferred array for blood typing may contain the following antibodies: anti-A, anti-B, anti-AB, anti-D (e.g., F8D8), anti-D (e.g., GAMA-401) and a suitable control (ligand only).

The following examples are provided so as to enable those of ordinary skill in the art to practice the method of the invention. These examples are not intended to limit the scope of what the inventors regards as their invention. Efforts have been made to ensure accuracy with respect to numbers used to characterize the measured conditions; however, some experimental errors and deviations may be present.



## EXAMPLE 1

Testing for the Presence or Absence of the A  
Bloodgroup Antigen on Erythrocytes using a Direct Assay

## 5 A. Preparation of Specific Anti-A Agarose Beads

Two ml of ImmunoPure® Immobilized Protein G was centrifuged for about two minutes at 1000 g and the supernate removed and discarded. Two ml of a pre-diluted sample of a mouse monoclonal anti-A (clone 46G5 mouse IgG3 isotype) was added. The monoclonal anti-A was diluted in 0.01M citrate buffered  
10 saline (pH 5.0) containing 7.5% Dextran (w/v). The contents were mixed gently and allowed to stand for 15 minutes. The agarose/antibody mixture was centrifuged for approximately two minutes at 1000g and the supernate removed and tested for residual anti-A activity in an agglutination assay.

15 That testing was performed in the following manner: a sample of the supernate was tested for its ability to agglutinate A positive erythrocytes. No agglutination was observed, indicating complete binding of the antibodies to the particles.

## 20 B. Preparation of the Reaction Tubes

Individual microwells of a microliter plate (Thermowell™ polycarbonate microtiter plate Costar) were used as reaction tubes. The specific anti-A agarose beads were thoroughly mixed and 70 $\mu$ l of the mixture was pipetted into each microwell. The agarose beads were allowed to settle to the bottom of the well  
25 before further testing (approximately 15 minutes).

## C. Test Procedure

Erythrocytes were obtained from either citrated blood samples or liquid nitrogen storage. The blood samples had been previously typed using FDA  
30 approved anti-A and anti-AB bloodgrouping reagents. A 0.8% suspension of the sample erythrocytes was prepared using a low ionic strength solution (Gamma

N-hance™, Gamma Biologicals, Inc.). Fifteen microliters (15 $\mu$ l) of the approximate 0.8% suspension of erythrocytes were pipetted into the reaction tube. The microtiter plate was then centrifuged for 15 seconds at 900-1000 g, 30 seconds at 500 g, and finally for 45 seconds at 900-1000 g in a Sero-fuge II (Clay-Adams), using a modified rotor so that the centrifugal force acts along the long axis of the reaction tube.

#### D. Results/Interpretation

If the erythrocytes possess the A bloodgroup antigen, they will adhere to the anti-A agarose beads at the top of the matrix (FIG. 1A). If the erythrocytes possess a weak expression of the A bloodgroup antigen, e.g., A<sub>3</sub> or A<sub>x</sub>, then some erythrocytes will adhere to the anti-A agarose beads at the top of the matrix and some will collect at the bottom of the wells (FIG. 1B). If the erythrocytes do not possess the A antigen, all the erythrocytes will collect at the bottom of the wells (FIG. 1C). The results observed in this experiment are shown below in Table 1.

TABLE 1  
ABO BLOODGROUP

	A <sub>1</sub>	A <sub>2</sub>	A <sub>1</sub> B	A <sub>2</sub> B	A <sub>3</sub>	A <sub>x</sub>	B	O
No. Tested	10	8	6	2	2	2	10	10
No. Positive	10	8	6	2	2	2	0	0

Testing erythrocytes for the presence or absence of the B bloodgroup antigen can be performed in a similar manner by substituting an anti-B for the anti-A reagent.

The source of the antibody can be either human polyclonal, human monoclonal or mouse monoclonal, etc. If the immunoglobulin class of the antibody is not IgG, then the antibody can be attached to a Protein G ligand by

means of a second antibody, e.g., anti-human IgM or anti-kappa or anti-lambda, provided that this secondary antibody is of the IgG class. Alternatively, if Protein A/G is bound to the gel particles, any antibody class to which Protein A/G binds may be used.

5

## EXAMPLE 2

### Testing for the Presence or Absence of the Rh D Antigen on Erythrocytes by Direct Assay

10

The erythrocytes were obtained from either citrated blood samples or liquid nitrogen storage. The erythrocytes had been previously typed using FDA approved anti-D bloodgrouping reagents. A 0.8% suspension of the erythrocytes was prepared in a low ionic strength solution (Gamma N-HANCE®, Gamma Biologicals, Inc.).

15

The preparation of the reaction tubes and test procedure are identical to those previously described in Example I except that an IgG1 anti-D (human monoclonal clone F8D8) was substituted for the anti-A antibodies.

20

#### A. Results/Interpretations

25

If the erythrocytes possess a normal expression of the rhesus D antigen, they will adhere to the anti-D agarose particles at the top of the solid phase matrix (FIG. 1A). If the erythrocytes possess a weak (quantitative) expression of the D antigen, e.g., D<sup>u</sup>, or a qualitative variant of the D antigen, e.g., DVI variant, then some of the erythrocytes will adhere to the anti-D agarose particles at the top of the matrix and some will collect at the bottom of the wells (FIG. 1B). If the erythrocytes do not possess the rhesus D antigen, all the erythrocytes will collect at the bottom of the wells (FIG. 1C). The results observed in this experiment are shown below in Table 2.

TABLE 2

Rh D

	Rh D Positive	Rh D Negative	DVI Variant	Weak D (D <sup>u</sup> )
No. Tested	10	6	3	3
No. Positive	10	0	3	2

Testing erythrocytes for the presence or absence of other Rh antigens, e.g., C, c, E or e, etc. can be performed in a similar manner by substituting that specific reagent for the anti-D antibodies. The source of the antibodies can be either human polyclonal (e.g., affinity purified), human monoclonal or mouse monoclonal, etc. If the immunoglobulin class of the antibody is not IgG, then the antibody can be attached to a Protein G ligand by means of a second antibody, e.g., anti-human IgM or anti-kappa or anti-lambda, provided that this secondary antibody is of the IgG class. Alternatively, if Protein A/G is bound to the gel particles, any antibody class to which Protein A/G binds may be used.

## EXAMPLE 3

Testing for the Presence or Absence of the  
Fy<sup>a</sup> Antigen on Erythrocytes by Indirect Assay

This test procedure was performed with non-affinity purified human polyclonal antisera.

## A. Preparation of Protein G Agarose Particles

Two ml of immobilized Protein G (ImmunoPure® G) was centrifuged for about two minutes at 900-1000g and the supernate removed and discarded. Two ml of a 0.01M citrate buffered saline (pH 5.0) solution containing 7.5% Dextran (w/v) was added to the Protein G and mixed thoroughly.

#### B. Preparation of the Reaction Tubes

Seventy microliters ( $70\mu\text{l}$ ) of the Protein G mixture from above was pipetted into microwells in a microtiter plate (Thermowell™ polycarbonate microtiter plate Costar). The agarose particles were allowed to settle to the base  
5 of the wells (approximately 15 minutes).

#### C. Test Procedure

Fifty microliters ( $50\mu\text{l}$ ) of a low ionic strength enhancement solution (Gamma N-HANCE®, Gamma Biologicals, Inc.) were added to each well. The  
10 erythrocytes were obtained either from citrated blood samples or from commercially available reagent erythrocyte panels (Gamma Biologicals, Inc.) and had been previously typed using FDA approved anti-Fy<sup>a</sup> bloodgrouping reagent. A 0.8% suspension of the erythrocytes was prepared using a low ionic strength solution. Fifteen microliters ( $15\mu\text{l}$ ) of the 0.8% suspension of the erythrocytes to  
15 be tested was added to each well followed by  $25\mu\text{l}$  of an anti-Fy<sup>a</sup> reagent test serum. The microtiter plate was then incubated for 10 minutes at 37°C.

The reaction tube was then centrifuged for 15 seconds at 900-1000g, 30 seconds at 500g, and 45 seconds at 900-1000g in a Sero-fuge II (Clay-Adams)  
20 using a modified rotor so that the centrifugal force acts directly along the axis of the reaction tube.

#### D. Results/Interpretation

If the erythrocytes possess the Fy<sup>a</sup> bloodgroup antigen, they will adhere to  
25 the Protein G agarose particles at the top of the matrix (FIG. 1A). If the erythrocytes do not possess the Fy<sup>a</sup> bloodgroup antigen, they will collect at the bottom of the wells (FIG. 1C). The results observed on this experiment are shown below in Table 3.

TABLE 3

Fy<sup>a</sup> Antigen Typing Duffy Bloodgroup Phenotype

	Fy(a+b-)	Fy(a+b+)	Fy(a-b+)
No. Tested	10	10	8
No. Positive	10	10	0

This test procedure could be used to type erythrocytes for any erythrocyte antigen in combination with a polyclonal or monoclonal antiserum containing specific IgG antibodies to the erythrocyte antigen.

## EXAMPLE 4

## Antibody Screening or Compatibility

## Testing by Indirect Assay

The erythrocytes used were from commercially available reagent erythrocyte panels (Gamma Biologicals, Inc.) that had been previously typed for specific antigens using FDA approved bloodgrouping reagents, where available, by the manufacturer. The sera samples were obtained from patient samples referred to the consultation laboratory of Gamma Biologicals, Inc. Each serum sample was tested with at least one erythrocyte suspension expressing the antigen against which the antibody was directed (positive control), and with at least one erythrocyte suspension lacking the antigen recognized by the antibody (negative control).

The preparation of the reaction tubes and test procedure are similar to those described in Example 3 except that either erythrocytes of a known antigenic composition (antibody screening) or donor erythrocytes (compatibility testing) are used in combination with a serum specimen from a patient or donor.

**A. Results/Interpretation**

If the patient's or donor's serum specimen contains an IgG antibody or antibodies directed against an antigen or antigens present on the erythrocytes used in the test procedure, the erythrocytes will adhere to the Protein G agarose particles at the top of the matrix (FIG. 1A).

In the case of a weak antigen-antibody reaction, some erythrocytes will adhere to the Protein G agarose particles at the top of the matrix and some will collect at the bottom of the reaction tube (FIG. 1B). If the serum specimen does not contain antibodies the antigens present on the erythrocytes, all the erythrocytes will collect at the bottom of the reaction tube (FIG. 1C). The results observed in this experiment are shown below in Table 4.

TABLE 4

Antibody Specificity	No. of Examples Tested	No. Positive
Anti-D	1	1
Anti-c	1	1
Anti-e	1	1
Anti-K	4	4
Anti-Fy <sup>a</sup>	3	3
Anti-Fy <sup>b</sup>	2	2
Anti-Jk <sup>a</sup>	2	2
Anti-Jk <sup>b</sup>	2	2
Anti-Xg <sup>a</sup>	2	2

## EXAMPLE 5

Testing for the Presence or Absence of the  
Rh D Antigen on Erythrocytes by Direct Assay  
Using an IgM Human Monoclonal Anti-D

5

## A. Preparation of Specific IgM Anti-D Agarose Beads

Protein G will not directly bind IgM antibodies. Thus, in order to produce a gel capable of binding an IgM antibody, a bridging antibody may first be attached to the Protein G agarose beads. In the following example, a mouse  
10 monoclonal antibody specific for human kappa light chains was utilized to attach an IgM kappa light chain human monoclonal anti-D to the agarose beads.

Two ml of Immunopure® Immobilized Protein G was centrifuged for approximately two minutes at 1000g and the supernate removed and discarded.  
15 Two ml of a pre-diluted sample of a mouse monoclonal anti-human kappa light chain (clone 5F3 mouse IgG1 isotype) was added. The monoclonal antibody containing culture supernate was diluted in 0.01M citrate buffered saline (pH 5.0). The contents were mixed gently and allowed to stand for 15 minutes. The agarose/antibody mixture was centrifuged for approximately two minutes at 1000g  
20 and the supernate removed and tested for residual anti-kappa activity in an agglutination assay.

That testing was performed in the following manner: a sample of pre-diluted anti-kappa not incubated with the Protein G agarose beads and a  
25 sample of diluted anti-kappa that had been pre-incubated with the Protein G agarose beads were tested in parallel for their ability to agglutinate erythrocytes previously coated with an IgG kappa light chain antibody. The preincubated sample gave a negative reaction, indicating that all the anti-kappa had become attached to the Protein G agarose beads.

30



Two ml of a sample of a human monoclonal anti-D (clone GAMA-401 IgM kappa light chain) was added to the anti-kappa/Protein G agarose pellet. The monoclonal antibody containing culture supernate was diluted in 0.01M citrate buffered saline (pH 5.0). The contents were mixed gently and allowed to stand  
5 for 15 minutes. The agarose/antibody mixture was centrifuged for approximately two minutes at 1000g and the supernate removed and tested for residual anti-D activity in an agglutination assay.

That testing was performed in the following manner: a sample of anti-D not  
10 incubated with the anti-kappa/Protein G agarose beads and a sample that had been pre-incubated with the anti-kappa/Protein G agarose beads were titrated in parallel and tested for their ability to agglutinate Rh D positive erythrocytes. The pre-incubated sample produced a significant reduction in titration value, indicating that anti-D had become attached to the matrix.

15 The anti-D/anti-kappa/Protein G agarose beads were washed three times with 0.85% NaCl in order to remove any unbound anti-D. After the final wash, two ml of a 0.01M citrate buffered saline solution (pH 5.0) containing 7.5% Dextran (w/v) were added to the beads and mixed thoroughly.

20 The erythrocytes were obtained from citrated blood samples which had been previously typed using FDA approved anti-D bloodgrouping reagents. A 0.8% suspension of erythrocytes was prepared in a low ionic strength solution (Gamma N-HANCE®, Gamma Biologicals, Inc.).

25 The preparation of the reaction tubes and test procedure are substantially identical to those previously described in Example 1.

A. Results/Interpretation

30 If the erythrocytes possess a normal expression of the rhesus D antigen, they will adhere to the anti-D agarose particles at the top of the solid phase matrix

(FIG. 1A). If the erythrocytes possess a weak (quantitative) expression of the D antigen, e.g., weak D, formerly called D<sup>u</sup>, then some will adhere to the anti-D agarose particles at the top of the matrix and some will collect at the bottom of the reaction tube (FIG. 1B). If the erythrocytes do not possess the rhesus D antigen,  
5 all the erythrocytes will collect at the bottom of the reaction tube (FIG. 1C). The that were observed in this list are shown by low in Table 5.

TABLE 5

	Rh D Positive	Rh D Negative	Weak D
10 No. Tested	8	5	2
No. Positive	8	0	2

Testing erythrocytes for the presence or absence of other bloodgroup  
antigens with IgM antibodies can be performed in a similar manner.

15

\* \* \*

Additional advantages and modifications will be readily apparent to those skilled in the art. The invention in its broader aspects is therefore not limited to  
20 the specific details or representative examples described. Accordingly, departures may be made from the detail without departing from the spirit or scope of the disclosed general inventive concept.

**CLAIMS:**

1. A method for detecting a bloodgroup antigen on erythrocytes, comprising the steps of:

5

obtaining a sample of erythrocytes to be tested;

10

adding said sample to a reaction tube, said tube having a lengthwise axis and containing a reaction medium, said reaction medium comprising a plurality of particles which have immunoglobulin-binding ligands selected from the group consisting of Protein A, Protein G, Protein A/G or a universal kappa light chain binding protein coupled to the surface of said particles and antibodies specific for said bloodgroup antigen coupled to said ligand on said particles;

15

centrifuging said reaction tube for a time sufficient to pellet to the bottom of said reaction tube erythrocytes that do not attach to said antibodies on said particles;

20

detecting attachment of said erythrocytes to said particles or lack thereof; and

correlating attachment with the presence of said bloodgroup antigen.

25

2. A method for detecting a bloodgroup antigen on erythrocytes, comprising the steps of:

30

obtaining a sample of erythrocytes to be tested;

5 adding said sample to a reaction tube, said tube having a lengthwise axis and containing a reaction medium, said reaction medium comprising a plurality of particles which have immunoglobulin-binding ligands selected from the group consisting of Protein A, Protein G, Protein A/G or a universal kappa light chain binding protein coupled to the surface of said particles, a bridging antibody coupled to said ligand and antibodies specific for said bloodgroup antigen coupled to said bridging antibody on said particles;

10 centrifuging said reaction tube for a time sufficient to pellet to the bottom of said reaction tube erythrocytes that do not attach to said antibodies coupled to said bridging antibodies that are coupled to said ligands on said particles;

15 detecting attachment of said erythrocytes to said particles or lack thereof, and

correlating attachment with the presence of said bloodgroup antigen.

20 3. The method of claims 1 or 2, wherein said centrifuging is performed using a centrifuge having a rotor adapted so that a centrifugal force generated by said centrifuging acts along said axis of said reaction tube.

25 4. The method of claim 3, wherein said centrifuging is performed for a time ranging from about 10-30 seconds in a range from about 800-1200g, for a time ranging from about 15-45 seconds in a range from about 300-600g, and for a time ranging from about 45-90 seconds in a range from about 800-1200g.

30

5. A method for detecting a blood cell antigen in blood typing, comprising the steps of:

obtaining a sample of erythrocytes to be tested;

5

incubating said erythrocytes in a reaction tube with a preparation of antibodies specific for a known blood cell antigen under conditions and for a time sufficient to bind said antibodies to said antigens if said antigens are present, said tube having a lengthwise axis and also containing a plurality of particles which have immunoglobulin-binding ligands selected from the group consisting of Protein A, Protein G, Protein A/G or a universal kappa light chain binding protein coupled to the surface of said particles;

10

centrifuging said reaction tube for a time sufficient to pellet to the bottom of said reaction tube erythrocytes that do not attach to said ligand on said particles; and

15

correlating attachment of said erythrocytes to said particles with the presence of said antigen on said erythrocytes.

20

6. A method for detecting blood serum antibodies specific for blood cell antigens, comprising the steps of:

25

obtaining a sample of erythrocytes having known antigens on the surface thereof;

obtaining a sample of blood serum to be tested for antibodies against said antigen;

30

-28-

- incubating said erythrocytes and said serum in a reaction tube under conditions and for a time sufficient for said antibodies to bind to said antigens on said erythrocytes if said antibodies are present, said tube having a lengthwise axis and also containing a plurality of particles having a ligand comprising an immunoglobulin-binding protein selected from the group consisting of Protein A, Protein G, Protein A/G or a universal kappa light chain binding protein coupled to the surface of said particles;
- centrifuging said reaction tube for a time sufficient to pellet to the bottom of said reaction tube erythrocytes that do not attach to said ligand on said particles; and
- correlating attachment of said erythrocytes to said particles with the presence of said antibodies.

7. The method of claims I and 5, wherein said antibodies are selected from the group consisting of anti-A, anti-B, anti-AB, anti-D, anti-C, anti-c, anti-E, anti-e, anti-K, anti-Fy<sup>a</sup>, anti-Fy<sup>b</sup>, anti-Jk<sup>a</sup>, anti-Jk<sup>b</sup>, anti-S or anti-s.

8. The method of claims 5 and 6, wherein said centrifuging is performed using a centrifuge having a rotor adapted so that a centrifugal force generated by said centrifuging acts along the axis of said reaction tube.

9. The method of claim 8, wherein said centrifuging is performed for a time ranging from about 10-30 seconds in a range from about 800-1200g, for a time ranging from about 15-45 seconds in a range from about 300-600g, and for a time ranging from about 45-90 seconds in a range from about 800-1200g.

10. A method for detecting blood serum antibodies specific for blood cell antigens, comprising the steps of:

- 5 obtaining a sample of erythrocytes having known antigens on the surface thereof;
- obtaining a sample of blood serum to be tested for antibodies against said antigens;
- 10 incubating for about 10 minutes at about 37°C said erythrocytes and said serum in a reaction tube, said tube having a lengthwise axis and containing a plurality of particles having Protein G coupled to the surface of each of said particles;
- 15 centrifuging said reaction tube for about 15 seconds at from about 900-1000g, then for about 30 seconds at from about 500g, and then for about 45 seconds at from about 900-1000g using a centrifuge adapted so that a centrifugal force generated by said centrifuge acts
- 20 along the axis of said reaction tube;
- detecting attachment of said erythrocytes to said Protein G on said particles or lack thereof; and
- 25 correlating attachment with the presence of said antibodies tested for.

11. The method of claims 1, 2, 5, 6 and 10, wherein said particles are of a material selected from the group consisting of agarose, Sepharose or latex.

12. The method of claims 1, 2, 5, 6 and 10, wherein said erythrocytes are prepared for testing in a low ionic strength solution creating an about 0.5-1.0% (w/v) suspension of said erythrocytes.

5

13. The method of claims 6 and 10, wherein said known antigens are selected from a group consisting of A, B, D, C, c, C<sup>w</sup>, E, e, K, Fy<sup>a</sup>, Fy<sup>b</sup>, Jk<sup>a</sup>, Jk<sup>b</sup>, S and s.

10

14. The method of claims 6 and 10, wherein said antibody tested for is selected from the group consisting of anti-A, anti-B, anti-AB, anti-D, anti-C, anti-c, anti-E, anti-e, anti-K, anti-Fy<sup>a</sup>, anti-Fy<sup>b</sup>, anti-Jk<sup>a</sup>, anti-Jk<sup>b</sup>, anti-S or anti-s.

15

15. The method of claim 5, wherein the antigen tested for is selected from the group consisting of A, B, D, C, c, C<sup>w</sup>, E, e, K, Fy<sup>a</sup>, Fy<sup>b</sup>, JK<sup>a</sup>, JK<sup>b</sup>, S and s.

20

16. An apparatus useful for detecting bloodgroup antigens and antibodies, comprising:

25

a plurality of reaction tubes being spaced apart and coupled together to form a single unit array of said reaction tubes, said array being adapted for use in a centrifuge;

each of said reaction tubes having:

30

an upper longitudinal region having a substantially rectangular cross-section of first preselected diameter, forming a reaction well adapted to receive reagents and erythrocytes;



- 31 -

5 a lower longitudinal region having a substantially circular cross-section of a second preselected diameter less than said first preselected diameter, forming a tube portion containing a column of immunoreactive particles having a ligand selected from the group consisting of Protein A, Protein G, Protein A/G or a universal kappa light chain binding protein coupled to the surface of said particles; and

10 an intermediate longitudinal region having a diameter varying between said first and said second preselected diameters, said intermediate longitudinal region providing fluid communication between said upper and lower longitudinal regions.

15 17. The apparatus of claim 16, wherein said single unit array of reaction tubes is in a linear configuration and further comprises a top adapted to sealably close each of said reaction tubes in said single unit array.

20 18. The apparatus of claim 17, wherein said single unit array has at least one substantially planar surface adapted for affixing identifying information.

25 19. The apparatus of claim 18, wherein said single unit array of reaction tubes is made of a biologically and chemically inert plastic.

20. The apparatus of claim 19, wherein said single unit array of reaction tubes is made of polymethylpentene.

Fig. 1a



Fig 1b



**Fig 1c**



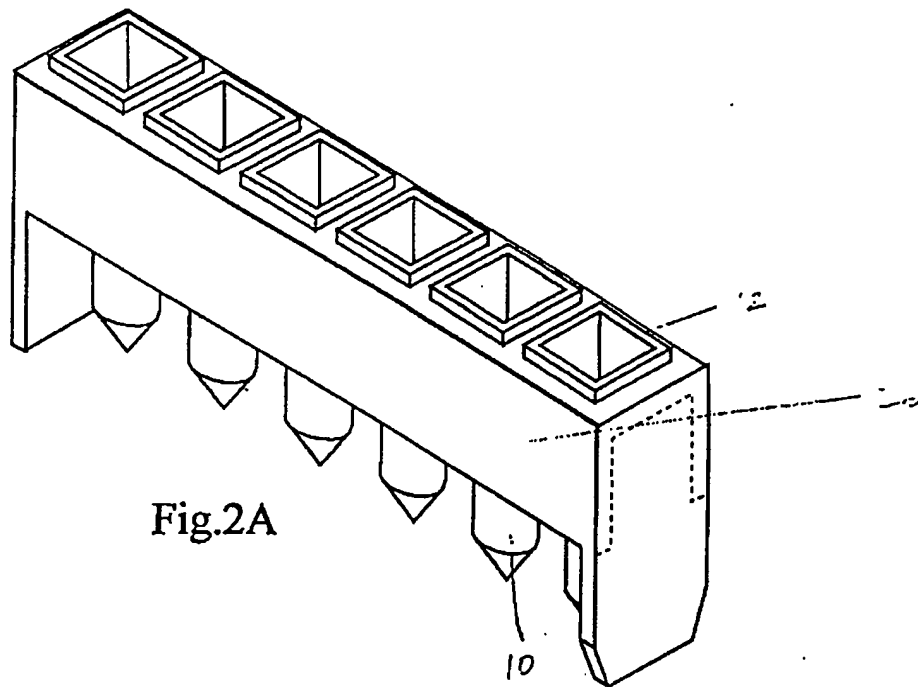


Fig. 2A

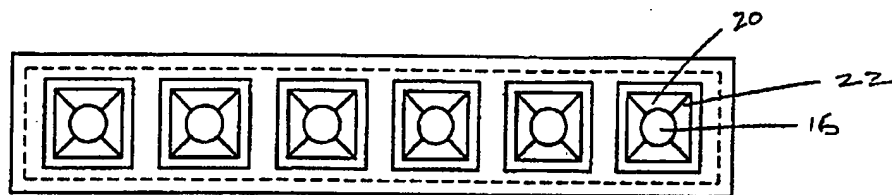


Fig. 2B

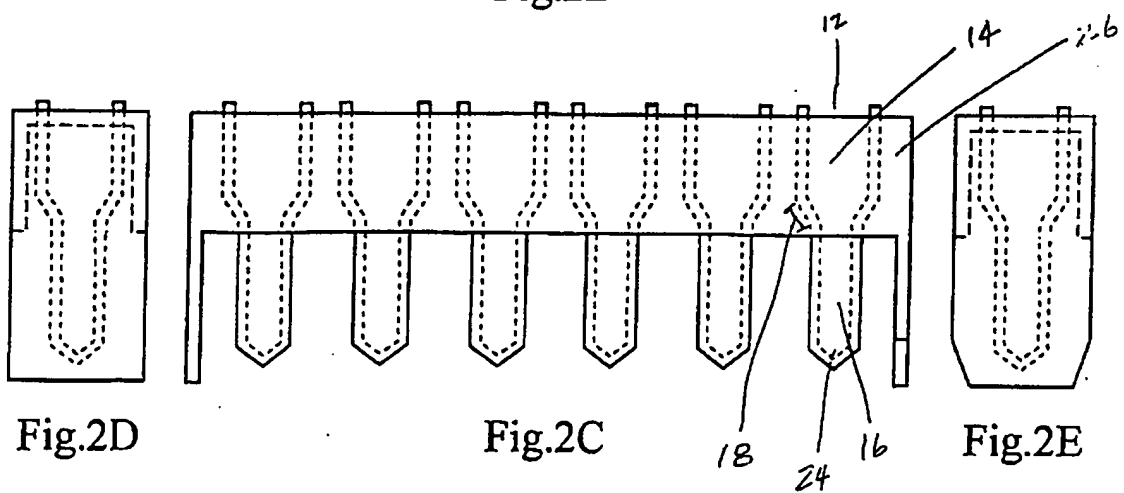


Fig. 2D

Fig. 2C

Fig. 2E

# INTERNATIONAL SEARCH REPORT

International Application No  
PCT/US 95/05982

A. CLASSIFICATION OF SUBJECT MATTER  
IPC 6 G01N33/80 G01N33/543

According to International Patent Classification (IPC) or to both national classification and IPC

## B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)  
IPC 6 G01N

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

## C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
A	EP,A,0 305 337 (STIFTUNG FUR DIAGNOSTISCHE FOR) 1 March 1989 ---	
A	EP,A,0 485 228 (ORTHO DIAGNOSTIC SYSTEMS INC) 13 May 1992 ---	
A	EP,A,0 512 896 (PASTEUR SANOFI DIAGNOSTICS) 11 November 1992 ---	
A	TRANSFUSION, vol. 30, 1990 pages 109-113, Y. LAPIERRE ET AL. 'The gel test: a new way to detect red cell antigen-antibody reactions.' cited in the application -----	

☐ Further documents are listed in the continuation of box C.

☒ Patent family members are listed in annex.

### \* Special categories of cited documents :

- \* "A" document defining the general state of the art which is not considered to be of particular relevance
- \* "E" earlier document but published on or after the international filing date
- \* "L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)
- \* "O" document referring to an oral disclosure, use, exhibition or other means
- \* "P" document published prior to the international filing date but later than the priority date claimed

- \* "T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
- \* "X" document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
- \* "Y" document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.
- \* "&" document member of the same patent family

Date of the actual completion of the international search

25 September 1995

Date of mailing of the international search report

10. 10. 95

Name and mailing address of the ISA  
European Patent Office, P.B. 5818-Patentlaan 2  
NL - 2280 HV Rijswijk  
Tel. (+ 31-70) 340-2040, Tx. 31 651 epo nl,  
Fax: (+ 31-70) 340-3016

Authorized officer

Cartagena y Abella,P

# INTERNATIONAL SEARCH REPORT

International Application No  
PCT/US 95/05982

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
EP-A-0305337	01-03-89	DE-A- 3875202	12-11-92
		JP-A- 1163662	27-06-89
		US-A- 5338689	16-08-94
-----			
EP-A-0485228	13-05-92	CA-A- 2055095	10-05-92
		GR-A- 91100453	08-10-92
		JP-A- 4285858	09-10-92
-----			
EP-A-0512896	11-11-92	FR-A- 2676123	06-11-92
		AT-T- 112392	15-10-94
		DE-D- 69200460	03-11-94
		DE-T- 69200460	11-05-95
		ES-T- 2061324	01-12-94
		JP-A- 5264543	12-10-93
		US-A- 5302512	12-04-94
-----			

**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ **BLACK BORDERS**
- ☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- ☐ **FADED TEXT OR DRAWING**
- ☐ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- ☐ **SKEWED/SLANTED IMAGES**
- ☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- ☐ **GRAY SCALE DOCUMENTS**
- ☒ **LINES OR MARKS ON ORIGINAL DOCUMENT**
- ☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- ☐ **OTHER:** \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**